



Medtronic


EOPA™ ELONGATED ONE-PIECE ARTERIAL CANNULA



Elongated for Easier Access During Minimally Invasive Procedures


EOPA™ ELONGATED ONE-PIECE ARTERIAL CANNULA

The Elongated One-Piece Arterial Cannula from Medtronic features a longer body for easier access during minimally invasive procedures.



Features	Benefits
<ul style="list-style-type: none"> ■ Vented plug ■ Introducer ■ Hemostasis cap ■ 3/8" vented or non-vented connectors 	<ul style="list-style-type: none"> ■ Maintains Hemostasis of the guidewire port ■ Maintains cannulae Hemostasis during insertion and priming ■ Flow-Guard™ style bloodless introduction ■ Vent facilitates de-airing after connection
<ul style="list-style-type: none"> ■ 13" length (33.02 cm) 	<ul style="list-style-type: none"> ■ Enables access to most areas of aortic arch through median sternotomy or peri sternotomy incisions
<ul style="list-style-type: none"> ■ Thin wall construction 	<ul style="list-style-type: none"> ■ Maximizes flow rates
<ul style="list-style-type: none"> ■ Wirewound body 	<ul style="list-style-type: none"> ■ Kink resistant, maximizes flexibility
<ul style="list-style-type: none"> ■ 18, 20, 22, 24 Fr. sizes (6.0 mm, 6.7 mm, 7.3 mm, 8.0 mm.) 	<ul style="list-style-type: none"> ■ Selection of sizes
<ul style="list-style-type: none"> ■ Adjustable radiopaque suture ring 	<ul style="list-style-type: none"> ■ Allows pre-setting of insertion depth and location for ligature sutures
<ul style="list-style-type: none"> ■ Depth marks 	<ul style="list-style-type: none"> ■ Indicate insertion depth at 1 cm intervals
<ul style="list-style-type: none"> ■ Flexible introducer with guidewire port 	<ul style="list-style-type: none"> ■ Facilitates insertion over 0.038" guidewire

Flow Performance

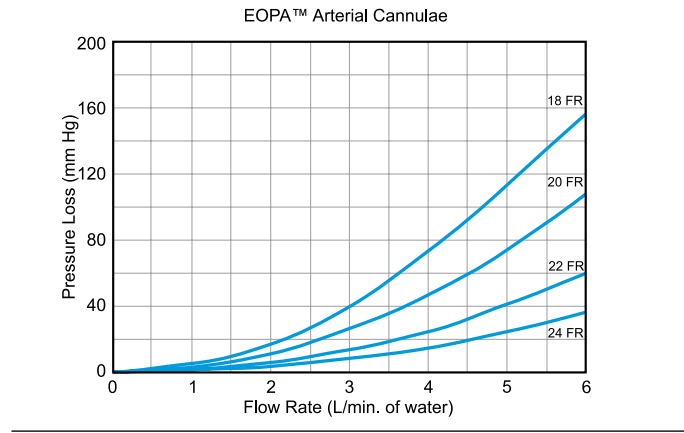

ORDERING INFORMATION

Catalog Code	Description	10/ctn.
77418	18 Fr. w/Vented 3/8" Connector	
77420	20 Fr. w/Vented 3/8" Connector	
77422	22 Fr. w/Vented 3/8" Connector	
77424	24 Fr. w/Vented 3/8" Connector	
77518	18 Fr. w/Non-Vented 3/8" Connector	
77520	20 Fr. w/Non-Vented 3/8" Connector	
77522	22 Fr. w/Non-Vented 3/8" Connector	
77524	24 Fr. w/Non-Vented 3/8" Connector	
77618	18 Fr. w/Dilator Introducer and Guidewire, 3/8" Vented Connector	
77620	20 Fr. w/Dilator Introducer and Guidewire, 3/8" Vented Connector	
77622	22 Fr. w/Dilator Introducer and Guidewire, 3/8" Vented Connector	
77624	24 Fr. w/Dilator Introducer and Guidewire, 3/8" Vented Connector	
77718	18 Fr. w/Dilator Introducer and Guidewire, 3/8" Non-Vented Connector	
77720	20 Fr. w/Dilator Introducer and Guidewire, 3/8" Non-Vented Connector	
77722	22 Fr. w/Dilator Introducer and Guidewire, 3/8" Non-Vented Connector	
77724	24 Fr. w/Dilator Introducer and Guidewire, 3/8" Non-Vented Connector	



World Headquarters
 Medtronic, Inc.
 710 Medtronic Parkway NE
 Minneapolis, MN 55432-5604
 USA
 Internet: www.medtronic.com
 Telephone: (763) 515-4000
 FAX: (763) 514-4879

Medtronic USA, Inc.
 Toll-free: 1-800-328-2518

Medtronic
Cardiac Surgical Products
 620 Watson, SW
 Grand Rapids, MI 49504
 USA
 Telephone: (616) 242-5200
 Toll-free: 1-800-253-1540
 FAX: (616) 242-5262

Europe
 Medtronic Europe S.A.
 Route du Molliau
 CH-1131 Tolochenaz
 Switzerland
 Internet: www.medtronic.co.uk
 Telephone: (41 21) 802-7000
 FAX: (41 21) 802-7900

Caution: Federal law (USA) restricts this device to sale by or on the order of a physician.

For a listing of indications, contraindications, precautions and warnings, please refer to the Directions For Use.

This device is for single use only.

For more information, contact your local Medtronic Cardiac Surgical Products Sales Representative or call Customer Service toll-free at 1-800-328-1357.